

THERAPYPOOLS


PRODUCT DESCRIPTION

The efficiency of hydrotherapy is widely recognised by medical professionals, and at Variopool too, we are familiar with this popular form of aquatherapy. The use of movable pool floors is a must to guarantee the best results, as they allow the various target groups to make optimal use of the therapy pool facilities, and they ensure the water temperature can be efficiently regulated.

The movable floor can be equipped with movable stairs, a lift for people with disabilities or a ramp to allow rehabilitation patients to access the pool unaided. Underwater treadmills can also be installed, so that patients do not need to turn back once they reach the wall, boosting the efficiency of their workout.

For more product information please visit our new website www.variomedic.nl

REFERENCES

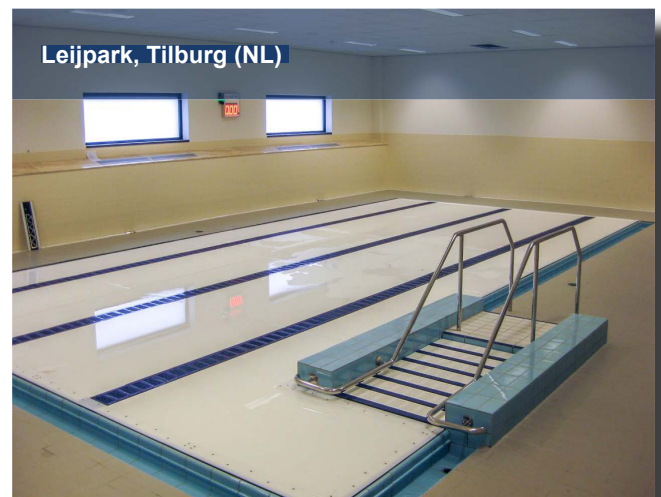

TECHNICAL SPECIFICATIONS

| | |
|----------------------------------|----------------------------------------------------------------------------------|
| Movable floor size | Project related |
| System | Buoyancy construction |
| Loadings | Min. 61,2 kg/m ² (EN13451-11) |
| Installation height | Max. 0.4mtr (50m ² max.) |
| Standard adjustment range | +0.05 to - 2.0mtr |
| Speed | 25 - 33 cm per minute |
| Material topdeck | Polyester, Polypropylene colour white RAL 9010 (standard)* |
| Material frame | Stainless steel 316L pickeld and passivated <i>* other colours on request</i> |

Below are some references where Variopool has designed and installed a movable floor.

Please contact us should you have questions about these or any other Variopool reference. You can reach us every day by calling +31 226 333 000 or sending an email to info@variomedic.nl

OPTIONS

Below are some references where Variopool has designed and installed a polypropylene movable pool floor. Please contact us should you have questions about these or any other Variomedic reference. You can reach us every day by calling +31 226 333 000 or sending an email to info@variomedic.nl

